


Baha[®] 6 Max Sound Processor



BAHA6 Max Sound Processor	Latest review date: May-2024	Version number: 1.6
Status: Approved	Next review: May-2026	Page 1 of 5

Components

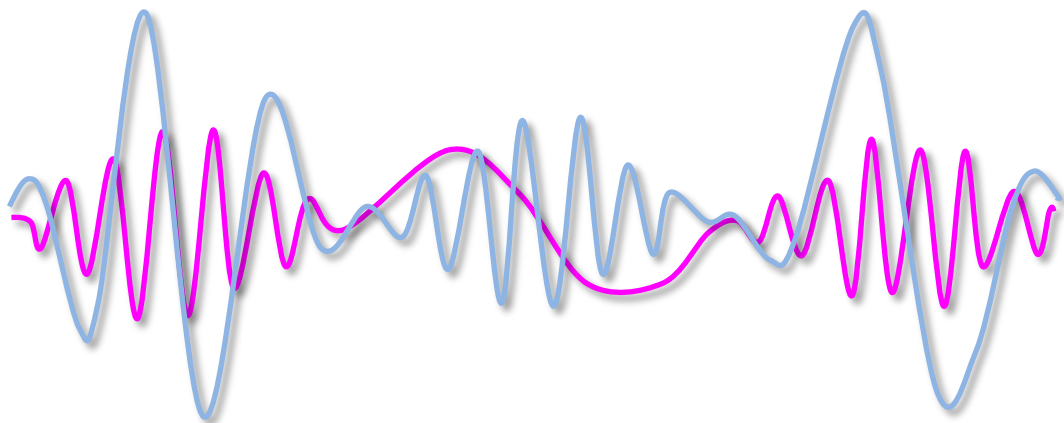

Tamper-proof battery door:

1. To unlock the device, carefully insert the tip of a pen in the small hole on the battery door and gently open the battery compartment. If the safety line is attached, one corner of the crocodile clip can also be used.



2. To lock the device, gently close the battery compartment until it is completely closed

When the battery is inserted, and the battery door is fully closed the Baha should produce several beeps before the microphone turns on. If the Baha is not positioned on the head when the microphone turns

BAHA6 Max Sound Processor	Latest review date: May-2024	Version number: 1.6
Status: Approved	Next review: May-2026	Page 2 of 5

on, expect to hear a buzzing / chirping noise (feedback) from it. A similar sound will be heard when the Baha is being worn and anything gets in close proximity to it.

Baha Safety line:



Note: The safety line should always be used when a Baha sound processor is worn by children.

Baha Soft band:



Features

1. Elasticated headband (free from natural rubber latex)
2. Universal slider
3. Built-in safety release
4. Snap connector
5. Baha Soft Wear Pad

Place the Baha Soft band on the wearer's head and adjust it to fit comfortably (it should be possible to insert one finger underneath the band). For best performance the snap connector should be placed flat behind the ear or at any other bony location of the skull i.e., forehead.

BAHA6 Max Sound Processor	Latest review date: May-2024	Version number: 1.6
Status: Approved	Next review: May-2026	Page 3 of 5

Attaching a Baha Soft Wear Pad:



You should change the Baha Soft Wear pad in the following cases:

- If the Baha Soft Wear pad appears worn-out or damaged
- If an accumulation of dirt or moisture on the Baha Soft Wear pad that cannot be wiped off
- If you experience a gradual deterioration in wearing comfort (if changing the Baha Soft Wear pad does not help, contact your hearing care professional)

Step-by-step

If necessary, remove the old pad from the snap connector. Ensure you remove any adhesive residue left behind.

1. Remove the cover sheet to expose the adhesive surface.
2. Attach the pad to the snap connector. Make sure it is centralized and covers the entire surface.
3. Remove the remaining cover sheet from the pad.

BAHA6 Max Sound Processor	Latest review date: May-2024	Version number: 1.6
Status: Approved	Next review: May-2026	Page 4 of 5

Attaching and removing:

To attach the sound processor onto the black disc, tilt it and gently snap it into place. To remove the sound processor place one finger underneath it and gently tilt until it releases.

Wireless Accessories are available and further information about these products and how they work can be found at:

<http://www.cochlear.com/wps/wcm/connect/uk/home/discover/baha-bone-conduction-implants/accessories-for-true-wireless-freedom>

<http://www.cochlear.com/wps/wcm/connect/uk/home/discover/cochlear-wireless-accessories/wireless-interactive-demo>

In the event the Baha develops a fault or you simply have any queries, please contact the Baha team directly cse.chimebaha@nhs.net.

If your personal details change, please inform us directly.

Address: Chime Social Enterprise, Audiology Department,
Royal Devon & Exeter Hospital (Wonford),
Barrack Road, Exeter, Devon, EX2 5DW

Tel: 01392 402223

E-Mail: cse.audiology@nhs.net

www.chimehealth.co.uk

Chime provides NHS Audiology Services for NHS Devon - Mid, East and Exeter areas. From 1st May 2011 the existing audiology department including staff and equipment transferred to Chime - the new Social Enterprise Company. Chime is located at the RD&E Foundation Trust (but is separately managed).

BAHA6 Max Sound Processor	Latest review date: May-2024	Version number: 1.6
Status: Approved	Next review: May-2026	Page 5 of 5